

VistaIntra DC – the smart choice



- Fast installation
- Perfectly matched to the image plate and sensor
- Constant radiation
- Variably adjustable mA and kV values
- DC generator with 0.4 mm focal spot
- Well thought-out options

Has your old X-ray emitter stopped working properly? The smart solution: VistalIntra DC



1. Sample old device



2. Old device removed



3. Simple installation of the adapter plate and electronics for the VistalIntra DC




4. The quick solution - VistalIntra DC

Fast & easy installation of adapter plates



2202-303-51
for Sirona
Heliodent DS

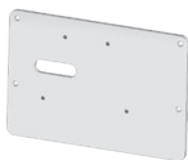

2202-303-52
for Sirona
Heliodent MD



2202-303-53
for Trophy
Irix CCX



2202-303-54
for Planmeca
Intra



2202-303-55
for Trophy Elity
(horizontal installation)



2202-303-50 Wall installation plates
for mounting on unstable walls

Well thought-out options



Cable adapter for connection of a 2-pin, potential-free push button
(included in the scope of delivery)



Remote Panel Kit for installation of the control panel outside the treatment room
(optional)



Image plate and film holder system Set Color

Fully equipped

Scope of delivery & options

Order no.	Product description
2202-01	VistaIntra DC with short extension arm Intraoral emitter with 200 mm tube and max. extension arm length of 1670 mm, Scissor-type arm, control panel with manual exposure switch. X-ray field collimation 3 x 4cm
2202-02	VistaIntra DC with medium extension arm Same as 2202-01, but with a max. extension arm length of 1894 mm
2202-03	VistaIntra DC with long extension arm Same as 2202-02, but with a max. extension arm length of 2194 mm
2202-101-50	X-ray field collimation 2 x 3cm Lead window 2 x 3cm
2202-101-51	X-ray field collimation 3 x 4cm Lead window 3 x 4cm (included in scope of delivery)
2202-102-50	Tube extension 300 mm Extends the tube from 200 to 300 mm
2202-302-50	Cover panel for VistaIntra DC Required for external attachment of the control panel
2130-981-50	Image plate and film holder system Set Color
2121-060-54	Intra / Extra Digital test body For acceptance and consistency checks
2202-303-50	Wall mounting plate Universal, for unstable walls
2202-303-51	Sirona Heliodent DS adapter plate For simple replacement of a Heliodent DS
2202-303-52	Sirona Heliodent MD adapter plate For simple replacement of a Heliodent MD
2202-303-53	Trophy IRIX adapter plate For simple replacement of a Trophy IRIX
2202-303-54	Planmeca Intra adapter plate For simple replacement of a Planmeca Intra
2202-303-55	Trophy Elitys adapter plate For simple replacement of a Trophy Elitys (horizontal installation only)

VistaSystem – the full range for digital diagnostics



VistaPano



VistaIntra



VistaScan Mini



VistaScan Combi



VistaCam

Count on outstanding image quality from a single source! With X-ray units, image plate scanners and intraoral camera systems from Dürr Dental, you can be confident that you will have the best possible diagnosis facilities offering maximum ease of use and the highest efficiency. Furthermore, our powerful imaging software perfectly links up all of the components of VistaSystem.